




Quick Aid
WHEN EVERY SECOND COUNTS

CELLTRION OTC RAPID ANTIGEN

Celltrion DiaTrust™

RAPID ANTIGEN

RESULTS IN ONLY 15-MINUTES

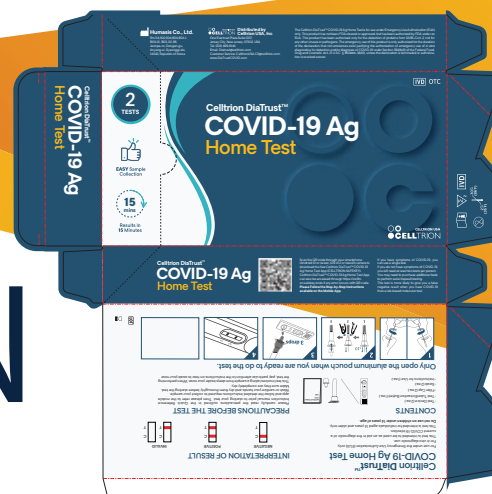

HIGHEST ACCURACY DETECTING COVID-19 & VARIANTS

- BEST STANDARD: FDA-EUA
- BEST RESULTS: 93.3% sensitivity / 99.3% specificity BEST
- PRICE: Please Inquire with your distributor
- BEST REVIEWS: 5-Stars for Safety & Ease of Use BEST
- AVAILABILITY: Stock on the ground now in USA!

The Other Tests...

ABBOTT Binax Now
84.6% sen / 98.5% spec
3-star Reviews
Low/Out of Stock

AccessBio CareStart
87.2% sen / 99.3% spec
4-star Reviews
Low/Out of Stock

QUIDEL QuickVue
78% sen / 97.5% spec
2-star Reviews
Low/Out of Stock



Quick Aid

WHEN EVERY SECOND COUNTS

CELLTRION OTC RAPID ANTIGEN

Celltrion DiaTrust™ COVID-19 Ag Home Test

INSTRUCTIONS FOR USE


For use under the Emergency Use Authorization (EUA) only. For *in vitro* diagnostic use. This test is intended to be used as an aid in the diagnosis of a current infection with the virus that causes COVID-19. This test is intended for individuals aged 14 years and older only. **Do not use on children under 14 years of age.**

STORAGE & STABILITY

Only open the aluminum pouch when you are ready to do the test. Immediately perform the test after opening the pouch. An unopened test device should be stored at **2 - 30°C (36 - 86°F)**. If the tests were refrigerated, keep them at room temperature for 30 minutes prior to use. It is stable until the expiration date marked on the label.

WHAT IS INCLUDED IN THIS BOX?

*The actual size of the test device may differ from the image.

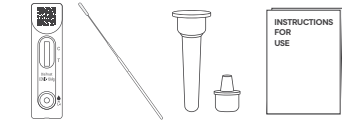

- ✓ Ensure all packaging is intact. Do not use the test if there is visible damage to the packaging or test pouch.

DOWNLOAD & OPEN APP

Scan the QR code through your smartphone (Android 10 or newer, iOS 14.2 or newer) camera to download the free Celltrion DiaTrust™ COVID-19 Ag Home Test App (CELLTRION SAFEKEY). Follow the instructions as described in the mobile app.



- ✓ Elderly population can acquire help from others to download & guide through the app.

Celltrion DiaTrust™ COVID-19 Ag Home Test App can also be accessed through <https://celltrion.safekey.tools> via a computer if any error occurs with QR code.


Please Follow the Step-by-Step Instructions Available on the Mobile App.

PRECAUTIONS BEFORE THE TEST

- ✓ Please carefully read the precautions outlined in the Instructions for Use manual prior to starting your test. Then please refer to the mobile app and follow the detailed instructions required to collect your sample. Failure to follow the instructions can result in inaccurate results.
- ✓ Wash or sanitize your hands and dry them thoroughly before starting the test. Make sure they are completely dry.
- ✓ This test involves taking a sample from deep inside your nose. When performing the test, pay particular attention to the instructions on how to swab your nose.
- ✓ Testing should be completed within 30-60 minutes of opening the test pouch.

TEST PROCEDURES


I. Swab Breakpoint and Holding Position

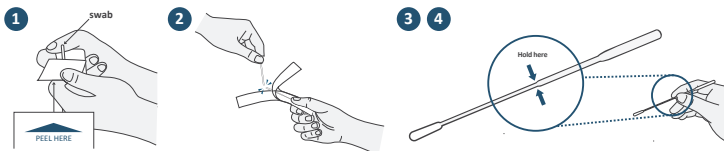

- Look for the "PEEL HERE" sign to peel open the swab package halfway. Make sure the soft tip is still covered with the packaging.
- Identify the breakpoint on the swab and break off the handle.
- Remove the swab from the package. Do not touch the soft tip or lay it down on any surfaces.
- You will see two notches on the handle. Make sure to hold the swab at the second notch, as pictured.

II. Collecting Your Nasal (mid-turbinate) Swab Sample

*Incorrect swabbing may lead to an inaccurate test result. This is particularly important if you do not have symptoms.

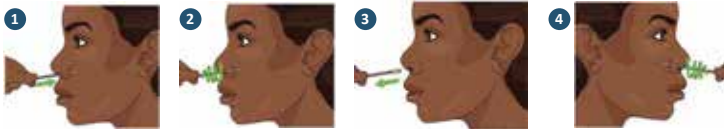


- Insert the entire soft end of the swab straight back into your nostril less than one inch (about 2 cm) or until resistance is felt.
- Slowly rotate the swab, gently rubbing it along the insides of your nasal passage several times.
- Gently remove the swab.
- Using the **same** swab, repeat this process in your other nostril with the same end of the swab.

Note: The swab included in the kit is designed for collection of samples from adults. Do not collect swabs from children under 14 years of age.

AFTER SAMPLE COLLECTION

- ✓ Put the tip of the swab into the test tube. Move the swab up and down at least 10 times to ensure sufficient sample extraction is extracted.
- ✓ Remove the swab while pressing against the sides of the tube to ensure maximum amount of liquid has been squeezed from the swab. **Note: False negative results can occur if the specimen is not properly mixed or too vigorously mixed.**
- ✓ Put the filter cap on the opening of the test tube and immediately dispense three drops of sample extract into the sample well of device. Click the "Completed" button and the 15 minute timer will start on the mobile application. **Note: Adding only one drop of solution or the entire vial may result in false negative results.**
- ✓ Read results at 15 minutes after applying the sample. Do not read results after 20 minutes. Click the "YES" or "NO" button on the application for the presence of red colored lines in the device window next to each of the two letters, C (Control) and T (Test). Follow the instructions based on your test result. **Note: False negative or false positive can occur if results are read before 15 minutes or after 20 minutes.**

Please refer to figure below for AFTER SAMPLE COLLECTION

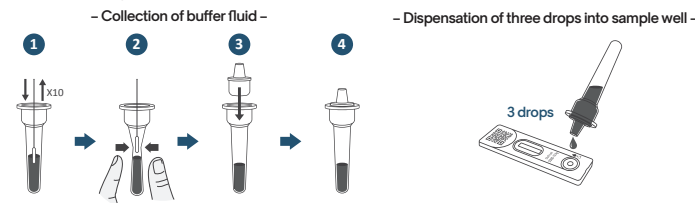

HOW TO READ THE RESULTS

A **Positive Result** indicates that viral antigens from COVID-19 were present in the specimen, and it is **very likely that you have COVID-19** and should self-isolate. It is important to be under the care of your healthcare provider. **Please make sure to compare your red colored line to the Line Level chart.**

Line Level	Result
0	NEGATIVE If no red colored line appears in the test line (T Line Level 0) and a red colored line is present on the control region (C Line Level 1 - 11), then the result is negative.
1-11	POSITIVE If a red colored line is visible in the test line (T Line Level 1 - 11) and control line (C Line Level 1 - 11), the results are positive.
12	INVALID If there is no red colored line in the control region (C Line Level 0), the result is invalid.

✓ If invalid results are obtained, please discard the device and repeat the testing using a **new** device.
✓ Dispose the test in general waste.

These are photos of actual positive results. Please note that the test line can show up faintly. This faint line still indicates a positive result.



- Negative results do not rule out COVID-19.
- In case of negative test result: Continue to follow all social distancing, recommendations and take protective measures. If suspicions of infection persist, repeat the test after 1 - 2 days and consult your doctor or local COVID-19 center.
- If your first test result is negative, you should test again in 24 to 48 hours.
- Note: A negative result is presumptive and additional testing with a molecular assay, may be needed.

The Celltrion DiaTrust™ COVID-19 Ag Home Test is for use under Emergency Use Authorization (EUA) only. This product has not been FDA cleared or approved, but has been authorized by FDA under an EUA. This product has been authorized only for the detection of proteins from SARS-CoV-2, not for any other viruses or pathogens. The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of *in vitro* diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated, or authorization is revoked sooner.



Quick Aid

WHEN EVERY SECOND COUNTS

CELLTRION OTC RAPID ANTIGEN

Celltrion DiaTrust™ COVID-19 Ag Home Test

INSTRUCTIONS FOR USE

INTENDED USE

Celltrion DiaTrust™ COVID-19 Ag Home Test is a lateral flow immunoassay intended for the qualitative detection of nucleocapsid protein and receptor binding domain (RBD) of the SARS-CoV-2 spike proteins in mid-turbinate swabs from the SARS-CoV-2 virus. This test is authorized for non-prescription home use with self-collected and adult-collected direct mid-turbinate swab samples from individuals aged 14 years or older with symptoms of COVID-19 within the first seven days of symptom onset.

This test is also authorized for non-prescription home use with self-collected and adult-collected mid-turbinate swab samples from individual 14 years or older with or without symptoms or other epidemiological reasons to suspect COVID-19 when tested twice over two or three days with at least 24 hours (and no more than 48 hours) between tests.

The Celltrion DiaTrust™ COVID-19 Ag Home Test does not differentiate between SARS-CoV and SARS-CoV-2.

Results are for the identification of SARS-CoV-2 nucleocapsid protein antigen and/or receptor binding domain (RBD). These antigens are generally detectable in mid-turbinate swabs during the acute phase of infection. Positive results indicate the presence of viral antigens, but clinical correlation with past medical history and other diagnostic information is necessary to determine infection status. Positive results do not rule out bacterial infection or co-infection with other viruses. The agent detected may not be the definite cause of disease. Individuals who test positive with the Celltrion DiaTrust™ COVID-19 Ag Home Test should self-isolate and seek follow up care with their physician or healthcare provider as additional testing and public health reporting may be necessary.

Negative results are presumptive, do not rule out SARS-CoV-2 infection and should not be used as the sole basis for treatment or patient management decisions, including infection control decisions. Negative results should be considered in the context of an individual's recent exposures, history and the presence of clinical signs and symptoms consistent with COVID-19 and confirmed with a molecular assay, if necessary, for patient management. For serial testing programs, additional confirmatory testing with a molecular test for negative results may be necessary, if there is a high likelihood of SARS-CoV-2 infection, such as an individual with a close contact with COVID-19 or with suspected exposure to COVID-19 or in communities with high prevalence of infection. Additional confirmatory testing with a molecular test for positive results may also be necessary, if there is a low likelihood of SARS-CoV-2 infection, such as in individuals without known exposures to COVID-19 or residing in communities with low prevalence of infection.

Individuals who test negative and continue to experience COVID-19 like symptoms of fever, cough and/or shortness of breath may still have SARS-CoV-2 infection and should seek follow up care with their physician or healthcare provider.

Individuals should provide all results obtained with this product to their healthcare provider for public health reporting. All healthcare providers will report all test results they receive from individuals who use the authorized product to relevant public health authorities in accordance with local, state, and federal requirements, using appropriate LOINC and SNOMED codes, as defined by the Laboratory In Vitro Diagnostics (LIVD) Test Code Mapping for SARS-CoV-2 Tests provided by CDC.

The Celltrion DiaTrust™ COVID-19 Ag Home Test is authorized for non-prescription self-use or an adult lay user testing another person 14 years or older. The Celltrion DiaTrust™ COVID-19 Ag Home Test is only for use under the Food and Drug Administration's Emergency Use Authorization.

WARNINGS & PRECAUTIONS

- Do not use this test for individuals under 14 years of age. The swab included in the kit is designed for collection of samples from adults and additional safety measures are needed for safe collection in children under 14 years of age.
- Testing for asymptomatic individuals should be performed at least twice over three days, with at least twenty-four hours and no more than 48 hours between tests. You

may need to purchase additional tests to perform this serial (repeat) testing.

- Do not use the test device beyond the expiration date.
- Keep sealed until usage, and once opened use immediately.
- Test samples immediately after collection.
- Do not use the test device if the test pouch is damaged or open.
- Do not re-use the device.
- This test is intended for diagnosis of coronavirus infection by detecting COVID-19 antigen, but should not be used as a sole criterion for the determination of SARS-CoV-2 infection. Other laboratory tests and clinical information (signs and symptoms) should be used and considered for diagnosis.
- Inadequate or inappropriate nasal swab sample collection may yield false test results.
- To obtain accurate results, the test must be performed as indicated in the application (Celltrion SafeKey) and/or instructions for use.
- Do not touch the swab head when handling the swab.
- Do not ingest extraction liquid.
- Keep out of reach of children.
- Avoid contact with skin and eyes.
- If contact with the body occurs, rinse with water. If irritation persists, seek medical advice.
- Discard Celltrion DiaTrust™ COVID-19 Ag Home Test in accordance with local, state and federal regulations or accreditation requirements.

FREQUENTLY ASKED QUESTIONS

WILL THIS TEST HURT?

No, the mid-turbinate nasal swab is not sharp, and it should not hurt. Sometimes the swab can feel slightly uncomfortable. If you feel pain, please stop the test and seek advice from your healthcare provider.

WHAT ARE THE KNOWN POTENTIAL RISKS AND BENEFITS OF THIS TEST?

- Potential risks include:
- Possible discomfort during sample collection.
 - Possible incorrect test results (see HOW TO READ THE RESULTS section).
- Potential benefits include:
- The results, along with other information, can help you and your healthcare provider make informed decisions about your care.
 - The results of this test may help limit the spread of COVID-19 to your family and others in your community.

For more information on EUAs go here: <https://www.fda.gov/emergency-preparedness-and-response/mcm-legal-regulatory-and-policy-framework/emergency-use-authorization>

WHAT IS SERIAL TESTING?

Serial testing is when a single person is tested for COVID-19 more than once. Because antigen tests are less sensitive than other COVID-19 tests and false results may occur, repeated testing may identify more individuals with COVID-19 infection than a single test. By repeating testing, it may be possible to more quickly identify cases of COVID-19 infection and reduce spread of infection. Additional testing with molecular COVID-19 test may be necessary, depending on your individual risk factors and test results. It is important that you work with your healthcare provider to help you understand the next steps you should take. Serial testing (i.e., testing every day or every other day) is more likely to detect COVID-19, especially when you do not have any symptoms.

WHAT IS THE DIFFERENCE BETWEEN AN ANTIGEN AND MOLECULAR TEST?

There are different kinds of tests for COVID-19. Molecular tests (also known as PCR tests) detect genetic material from the virus. Antigen tests, such as the Celltrion DiaTrust™ COVID-19 Ag Home Test detect proteins from the virus that causes COVID-19. Antigen tests are very specific for the COVID-19 virus but are not as sensitive as molecular tests. This means that a positive result is highly accurate, but a negative result does not rule out infection. If your test result is negative, you should discuss with your healthcare provider

whether an additional molecular test is necessary and if you should continue isolating at home. There is a higher chance of false negative results with antigen tests than with laboratory-based molecular tests. This means that there is a higher chance this test will give you a negative result when you have COVID-19.

HOW ACCURATE IS THIS TEST?

The clinical evaluation of the Celltrion DiaTrust™ COVID-19 Ag Home Test was evaluated by testing a total of 492 prospectively collected direct mid-turbinate nasal swab samples, consisted of 45 positive and 447 negative samples from suspected COVID-19 patients in United States that were within seven days of symptom onset or asymptomatic, aged 14 years and older. The Celltrion DiaTrust™ COVID-19 Ag Home Test was compared to an FDA authorized molecular SARS-CoV-2 test. The Celltrion DiaTrust™ COVID-19 Ag Home Test correctly identified 86.7% of positive specimens and 99.8% of negative specimens in that clinical study.

WHAT IF YOU TEST POSITIVE?

A positive test result means that antigens from COVID-19 were detected and it is very likely you currently have COVID-19 disease. There is a very small chance that this test can give a positive result that is wrong (a false positive result). If you test positive you should self-isolate at home per CDC recommendations to stop spreading the virus to others. Please consult the CDC recommendations regarding self-isolation at www.cdc.gov/coronavirus. Seek follow-up care with your healthcare provider immediately. Your healthcare provider will work with you to determine how best to care for you based on your test result(s) along with your medical history, and your symptoms.

WHAT IF YOU TEST NEGATIVE?

A negative test result indicates no antigens for COVID-19 were detected. It is possible for this test to give a negative result that is incorrect (false negative) in some people with COVID-19 and negative results are presumptive and may need to be confirmed with a molecular test. This means that you could possibly still have COVID-19 even though the test is negative. If you receive a negative result, you should test again in 24-48 hours. If you test negative and continue to experience symptoms of fever, cough and/or shortness of breath you should seek follow up care with your healthcare provider immediately. Your healthcare provider may suggest you need another test to determine if you have contracted the virus causing COVID-19. If you are concerned about your COVID-19 infection status after testing or think you may need follow up testing, please contact your healthcare provider.

HAZARDOUS INGREDIENT FOR REAGENT

Chemical Name (CAS)	GHS Code for each ingredient	Conc.
Sodium Azide (26628-22-8)	Acute Tox.2 (oral), H300 Acute Tox.1 (dermal), H310	0.09%

The extraction buffer solution in the extraction buffer tube contains a hazardous ingredient as shown in above table. If the extraction buffer solution contacts the skin or eye, immediately wash with plenty of running water. In case the irritation persists, please seek medical advice at: <https://www.poison.org/contact-us-or-1-800-222-1222>.

If you have any questions, please contact Humasis Co., Ltd. (via email: info@humasis.com, via phone: +82-31-9085-6284) or Celltrion USA, Inc. (via email: celltrionusa.CS@celltrion.com, or via phone: (201) 499-1844)

Humasis Co., Ltd.

Rm. 114, 502, 504, 604, 604-1, B03-1,
B03-2, 88, Jeonpa-ro, Dongan-gu,
Anyang-si, Gyeonggi-do, 14042,
Republic of Korea
Fax: (201) 603-6767

Distributed by Celltrion USA, Inc

One Evertrust Plaza Suite 1207,
Jersey City, New Jersey, 07302, USA
Tel: (201) 499-1844
Email: Diatrust@celltrion.com
www.DiaTrustCOVID.com



Quick Aid
WHEN EVERY SECOND COUNTS

CELLTRION OTC RAPID ANTIGEN



October 21, 2021

Ally Danta
Senior Associate
Parexel International
Representing: Celltrion USA, Inc.
2520 Meridian Parkway, Suite 200
Durham, NC 27713

Device: Celltrion DiaTrust COVID-19 Ag Home Test

EUA Number: EUA210501

Company: Celltrion USA, Inc.

Indication: Non-prescription home use for the qualitative detection of nucleocapsid protein antigen and receptor binding domain (RBD) of the SARS-CoV-2 spike proteins from SARS-CoV-2 with:

Self-collected and adult-collected direct mid-turbinate swab samples from individuals aged 14 years or older with symptoms of COVID-19 within the first seven days of symptom onset.

Self-collected and adult-collected mid-turbinate swab samples from individuals aged 14 years or older, with or without symptoms or other epidemiological reasons to suspect COVID-19 when tested twice over three days with at least 24 hours (and no more than 48 hours) between tests.

Dear Ms. Danta:

This letter is in response to your¹ request that the Food and Drug Administration (FDA) issue an Emergency Use Authorization (EUA) for emergency use of your product,² pursuant to Section 564 of the Federal Food, Drug, and Cosmetic Act (the Act) (21 U.S.C. §360bbb-3).

On February 4, 2020, pursuant to Section 564(b)(1)(C) of the Act, the Secretary of the Department of Health and Human Services (HHS) determined that there is a public health emergency that has a significant potential to affect national security or the health and security of United States citizens living abroad, and that involves the virus that causes COVID-19. Pursuant to Section 564 of the Act, and on the basis of such determination, the Secretary of HHS then declared that circumstances exist justifying the authorization of emergency use of in

¹ For ease of reference, this letter will use the term “you” and related terms to refer to Celltrion USA, Inc.

² For ease of reference, this letter will use the term “your product” to refer to the Celltrion DiaTrust COVID-19 Ag Home Test, manufactured by HUMASIS Co., Ltd., used for the indication identified above.



Quick Aid
WHEN EVERY SECOND COUNTS

CELLTRION OTC RAPID ANTIGEN

Page 2 – Ally Danta, Celltrion USA, Inc.

vitro diagnostics for detection and/or diagnosis of the virus that causes COVID-19 subject to the terms of any authorization issued under Section 564(a) of the Act.³

FDA considered the totality of scientific information available in authorizing the emergency use of your product for the indication above. A summary of the performance information FDA relied upon is included in the “Celltrion DiaTrust COVID-19 Ag Home Test Healthcare Provider Instructions for Use” identified below.

Having concluded that the criteria for issuance of this authorization under Section 564(c) of the Act are met, I am authorizing the emergency use of your product, described in the Scope of Authorization of this letter (Section II), subject to the terms of this authorization.

I. Criteria for Issuance of Authorization

I have concluded that the emergency use of your product meets the criteria for issuance of an authorization under Section 564(c) of the Act, because I have concluded that:

1. The SARS-CoV-2 can cause a serious or life-threatening disease or condition, including severe respiratory illness, to humans infected by this virus;
2. Based on the totality of scientific evidence available to FDA, it is reasonable to believe that your product may be effective in diagnosing COVID-19, and that the known and potential benefits of your product when used for diagnosing COVID-19, outweigh the known and potential risks of your product; and
3. There is no adequate, approved, and available alternative to the emergency use of your product.⁴

II. Scope of Authorization

I have concluded, pursuant to Section 564(d)(1) of the Act, that the scope of this authorization is limited to the indication above.

Authorized Product Details

Your product is a lateral flow immunoassay intended for the qualitative detection of nucleocapsid protein antigen and receptor binding domain (RBD) of the SARS-CoV-2 spike proteins from the SARS-CoV-2. This test is authorized for non-prescription home use with self-collected and adult-collected direct mid-turbinate swab samples from individuals aged 14 years or older with symptoms of COVID-19 within the first seven days of symptom onset.

This test is also authorized for non-prescription home use with self-collected mid-turbinate swab samples from individuals aged 14 years or older with or without symptoms or other

³ U.S. Department of Health and Human Services, *Determination of a Public Health Emergency and Declaration that Circumstances Exist Justifying Authorizations Pursuant to Section 564(b) of the Federal Food, Drug, and Cosmetic Act*, 21 U.S.C. § 360bbb-3. 85 FR 7316 (February 7, 2020).

⁴ No other criteria of issuance have been prescribed by regulation under Section 564(c)(4) of the Act.

CELLTRION OTC RAPID ANTIGEN

Page 3 – Ally Danta, Celltrion USA, Inc.

epidemiological reasons to suspect COVID-19 when tested twice over three days with at least 24 hours (and no more than 48 hours) between tests.

The SARS-CoV-2 nucleocapsid protein antigen and/or receptor binding domain (RBD) are generally detectable in mid-turbinate nasal swabs during the acute phase of infection. Positive results indicate the presence of viral antigens, but clinical correlation with past medical history and other diagnostic information is necessary to determine infection status. Positive results do not rule out bacterial infection or co-infection with other viruses. The agent detected may not be the definite cause of disease. Individuals who test positive with your product should self-isolate and seek follow up care with their physician or healthcare provider as additional testing and public health reporting may be necessary.

Negative results are presumptive, do not rule out SARS-CoV-2 infection and should not be used as the sole basis for treatment or patient management decisions, including infection control decisions. Negative results should be considered in the context of an individual's recent exposures, history and the presence of clinical signs and symptoms consistent with COVID-19 and confirmed with a molecular assay, if necessary, for patient management.

For serial testing programs, additional confirmatory testing with a molecular test for negative results may be necessary, if there is a high likelihood of SARS-CoV-2 infection, such as an individual with a close contact with COVID-19 or with suspected exposure to COVID-19 or in communities with high prevalence of infection. Additional confirmatory testing with a molecular test for positive results may also be necessary, if there is a low likelihood of SARS-CoV-2 infection, such as in individuals without known exposures to COVID-19 or residing in communities with low prevalence of infection.

Individuals who test negative and continue to experience COVID-19 like symptoms of fever, cough and/or shortness of breath may still have SARS-CoV-2 infection and should seek follow up care with their physician or healthcare provider.

Individuals should provide all results obtained with this product to their healthcare provider for public health reporting. All healthcare providers will report all test results they receive from individuals who use the authorized product to relevant public health authorities in accordance with local, state, and federal requirements using appropriate LOINC and SNOMED codes, as defined by the Laboratory In Vitro Diagnostics (LIVD) Test Code Mapping for SARS-CoV-2 Tests provided by the Centers for Disease Control and Prevention (CDC).

Your product is performed using mid-turbinate swab samples from individuals aged 14 years or older. The individual using your product is instructed to download, open and follow the step-by-step mobile- or Web-application based instructions on the Celltrion DiaTrust COVID-19 Ag Home Test Application (Celltrion SafeKey) onto a computer or compatible smartphone.⁵ When using your product the individual opens the aluminum pouches containing (1) the test device and (2) the test tube and filter cap. The QR code on the test device is then scanned with the

⁵ Compatible smart phone includes Apple iPhone running Operation System (iOS) 14.2 or later versions of the iOS, and Android Phones running Android 10 or later versions. Additional smart phone models as may be requested, and for which you receive appropriate authorization, in accordance with Condition Q. below.

CELLTRION OTC RAPID ANTIGEN

Page 4 – Ally Danta, Celltrion USA, Inc.

smartphone and the requested information completed. The swab is then removed from its packaging and the individual collects an mid-turbinate swab sample by inserting the swab into the nostril and slowly rotating the swab, gently rubbing it along the insides of the nasal passage several times, before repeating the process in the second nostril. The swab is then immediately inserted into the test tube and moved up and down at least 10 times. The swab is then removed while pressing against the sides of the tube and the tube capped with the filter cap. Three drops of the solution are applied into the sample well of the test device. The individual then clicks the “completed” button in the App that activates a 15 minute timer. When the extracted swab specimen is dispensed into to the sample well, the specimen migrates towards the conjugate pad, which contains conjugated antibodies with colloidal gold directed against the SARS-CoV-2 antigen. When the sample contains SARS-CoV-2 antigens, an antigen-antibody-conjugate complex is formed. The sample-conjugate complex then passes over the membrane until it reaches the capture zone (test line). Here, the complex is bound to immobilized antibodies and form visible colored band in the test line. The sample then migrates across the membrane along the strip until it reaches the control line where excess conjugate binds and produces a second visible line on the membrane. This control line indicates that the sample has migrated across the membrane as intended and indicates that the test was correctly performed. Test results are interpreted visually after 15 minutes based on the presence or absence of visually detectable colored lines at the control line (C) and/or test line (T). Upon completion of the test and result interpretation the user should share their results with their healthcare provider.

The Celltrion DiaTrust COVID-19 Ag Home Test includes the following materials or other authorized materials: Test Cassette with test strip, Extraction Buffer (provided in a sealed test tube), Filter cap, Swab and “Celltrion DiaTrust COVID-19 Ag Home Test Instructions for Use.”

Your product includes an internal control test line (“C”) that must generate the expected result for a test to be considered valid, as outlined in the “Celltrion DiaTrust COVID-19 Ag Home Test Instructions for Use” and the “Celltrion DiaTrust COVID-19 Ag Home Test Healthcare Provider Instructions for Use.”

The labeling entitled “Celltrion DiaTrust COVID-19 Ag Home Test Healthcare Provider Instructions for Use,” the “Celltrion DiaTrust COVID-19 Ag Home Test Instructions for Use,” and the “Celltrion DiaTrust COVID-19 Ag Home Test” box labels (available at <https://www.fda.gov/medical-devices/coronavirus-disease-2019-covid-19-emergency-use-authorizations-medical-devices/in-vitro-diagnostics-euas>), the “Celltrion DiaTrust COVID-19 Ag Home Test Application (Celltrion SafeKey)” software application and the following fact sheet pertaining to the emergency use, is required to be made available as set forth in the Conditions of Authorization (Section IV), and are collectively referred to as “authorized labeling”:

- Fact Sheet for Healthcare Professionals⁶: Celltrion USA, Inc. - Celltrion DiaTrust COVID-19 Ag Home Test

The above described product, when accompanied by the authorized labeling as set forth in the

⁶ Note that the information typically found in a Fact Sheet for Individuals is contained in the authorized “Celltrion DiaTrust COVID-19 Ag Home Test Instructions for Use ” that will be available to end users as set forth in the Conditions of Authorization (Section IV).

CELLTRION OTC RAPID ANTIGEN

Page 5 – Ally Danta, Celltrion USA, Inc.

Conditions of Authorization (Section IV) is authorized to be distributed and used under this EUA, despite the fact that it does not meet certain requirements otherwise required by applicable federal law.

I have concluded, pursuant to Section 564(d)(2) of the Act, that it is reasonable to believe that the known and potential benefits of your product, when used consistent with the Scope of Authorization of this letter (Section II), outweigh the known and potential risks of your product.

I have concluded, pursuant to Section 564(d)(3) of the Act, based on the totality of scientific evidence available to FDA, that it is reasonable to believe that your product may be effective in diagnosing COVID-19, when used consistent with the Scope of Authorization of this letter (Section II), pursuant to Section 564(c)(2)(A) of the Act.

FDA has reviewed the scientific information available to FDA, including the information supporting the conclusions described in Section I above, and concludes that your product (as described in the Scope of Authorization of this letter (Section II)) meets the criteria set forth in Section 564(c) of the Act concerning safety and potential effectiveness.

The emergency use of your product under this EUA must be consistent with, and may not exceed, the terms of this letter, including the Scope of Authorization (Section II) and the Conditions of Authorization (Section IV). Subject to the terms of this EUA and under the circumstances set forth in the Secretary of HHS's determination under Section 564(b)(1)(C) of the Act described above and the Secretary of HHS's corresponding declaration under Section 564(b)(1) of the Act, your product is authorized for the indication above.

III. Waiver of Certain Requirements

I am waiving the following requirements for your product during the duration of this EUA:

- Current good manufacturing practice requirements, including the quality system requirements under 21 CFR Part 820 with respect to the design, manufacture, packaging, labeling, storage, and distribution of your product, but excluding Subpart H (Acceptance Activities, 21 CFR 820.80 and 21 CFR 820.86), Subpart I (Nonconforming Product, 21 CFR 820.90), and Subpart O (Statistical Techniques, 21 CFR 820.250).

IV. Conditions of Authorization

Pursuant to Section 564(e) of the Act, I am establishing the following conditions on this authorization:

Celltrion USA, Inc. (You), HUMASIS Co., Ltd. and Authorized Distributor(s)⁷

⁷ "Authorized Distributor(s)" are identified by you, Celltrion USA, Inc., in your EUA submission as an entity allowed to distribute your product.



Quick Aid
WHEN EVERY SECOND COUNTS

CELLTRION OTC RAPID ANTIGEN

Page 6 – Ally Danta, Celltrion USA, Inc.

- A. Your product must comply with the following labeling requirements: the intended use statement in 21 CFR 809.10(a)(2), (b)(2); adequate directions for use in 21 U.S.C. 352(f) and 21 CFR 809.10(b)(5), (7), and (8); appropriate limitations on the use of the device including information required under 21 CFR 809.10(a)(4); and any available information regarding performance of the device, including requirements under 21 CFR 809.10(b)(12).
- B. You, HUMASIS Co., Ltd. and authorized distributor(s) must maintain records of customer complaint files and report to FDA any significant complaints about usability or deviations from the established performance characteristics of which you and authorized distributor(s) become aware.
- C. You, HUMASIS Co., Ltd. and authorized distributor(s) must collect information on the performance of your product and have a process in place to track adverse events, including any occurrence of false positive or false negative results and significant deviations from the established performance characteristics of the product of which you become aware and report any such events to FDA in accordance with 21 CFR Part 803. Serious adverse events, especially unexpected biosafety concerns, should immediately be reported to the Division of Microbiology (DMD)/Office of Health Technology 7 (OHT7)-Office of In Vitro Diagnostics and Radiological Health (OIR)/Office of Product Evaluation and Quality (OPEQ)/Center for Devices and Radiological Health (CDRH) (via email: CDRH-EUAREporting@fda.hhs.gov).
- D. You, HUMASIS Co., Ltd. and authorized distributor(s) using your product must ensure that any records associated with this EUA are maintained until otherwise notified by FDA. Such records will be made available to FDA for inspection upon request.

Celltrion USA, Inc. (You) and HUMASIS Co., Ltd.

- E. You and HUMASIS Co., Ltd. must comply with the following requirements pursuant to FDA regulations: 21 CFR 820 Subpart H (Acceptance Activities, 21 CFR 820.80 and 21 CFR 820.86), Subpart I (Nonconforming Product, 21 CFR 820.90), and Subpart O (Statistical Techniques, 21 CFR 820.250).
- F. You and HUMASIS Co., Ltd. must have lot release procedures and the lot release procedures, including the study design and statistical power, must ensure that the tests released for distribution have the clinical and analytical performance claimed in the authorized labeling.
- G. If requested by FDA, you and HUMASIS Co., Ltd. must submit lot release procedures to FDA, including sampling protocols, testing protocols, and acceptance criteria, that you use to release lots of your product for distribution in the U.S. If such lot release procedures are requested by FDA, you must provide it within 48 hours of the request.
- H. You and HUMASIS Co., Ltd. must evaluate the impact of SARS-CoV-2 viral mutations on your product's performance. Such evaluations must occur on an ongoing basis and



Quick Aid
WHEN EVERY SECOND COUNTS

CELLTRION OTC RAPID ANTIGEN

Page 7 – Ally Danta, Celltrion USA, Inc.

must include any additional data analysis that is requested by FDA in response to any performance concerns you or FDA identify during routine evaluation. Additionally, if requested by FDA, you must submit records of these evaluations for FDA review within 48 hours of the request. If your evaluation identifies viral mutations that affect the stated expected performance of your device, you must notify FDA immediately (via email: CDRH-EUA-Reporting@fda.hhs.gov).

- I. If requested by FDA, you and HUMASIS Co., Ltd. must update your labeling within 7 calendar days to include any additional labeling risk mitigations identified by FDA, such as those related to the impact of viral mutations on test performance. Such updates will be made in consultation with, and require concurrence of, DMD/OHT7-OIR/OPEQ/CDRH.

Celltrion USA, Inc. (You) and Authorized Distributor(s)

- J. You and authorized distributor(s) must make available the “Celltrion DiaTrust COVID-19 Ag Home Test Instructions for Use” for use in the shipped kit using the “Celltrion DiaTrust COVID-19 Ag Home Test” box labels and make these two documents electronically available on your website.
- K. You and authorized distributor(s) must inform relevant public health authorities of this EUA, including the terms and conditions herein, and any updates made to your product and/or the authorized labeling.
- L. Through a process of inventory control, you and authorized distributor(s) must maintain records of the locations (e.g., pharmacies, doctor’s offices, etc.) to which your product is distributed and the number of tests distributed to each location.
- M. You and authorized distributor(s) are authorized to make available additional information relating to the emergency use of your product that is consistent with, and does not exceed, the terms of this letter of authorization.

Celltrion USA, Inc. (You)

- N. You must notify FDA of any authorized distributor(s) of your product, including the name, address, and phone number of any authorized distributor(s).
- O. You must provide authorized distributor(s) with a copy of this EUA and communicate to authorized distributor(s) any subsequent revisions that might be made to this EUA and its authorized accompanying materials, including the authorized labeling.
- P. You must make the authorized “Celltrion DiaTrust COVID-19 Ag Home Test Healthcare Provider Instructions for Use” and the “Fact Sheet for Healthcare Professionals” electronically available on your website. Additionally, you must provide the opportunity to request a copy of the “Celltrion DiaTrust COVID-19 Ag Home Test Healthcare Provider Instructions for Use” and “Fact Sheet for Healthcare Professionals” in paper



CELLTRION OTC RAPID ANTIGEN

Page 8 – Ally Danta, Celltrion USA, Inc.

form, and after such request, promptly provide the requested labeling at no additional cost.

- Q. You may request changes to this EUA for your product, including to the Scope of Authorization (Section II in this letter) or to the authorized labeling, including requests to make available additional authorized labeling specific to an authorized distributor. Such additional labeling may use another name for the product but otherwise must be consistent with the authorized labeling, and shall not exceed the terms of authorization of this letter. Any request for changes to this EUA should be submitted to DMD/OHT7-OIR/OPEQ/CDRH and require appropriate authorization from FDA prior to implementation.
- R. You must evaluate the analytical limit of detection and assess traceability⁸ of your product with any FDA-recommended reference material(s). After submission to and concurrence with the data by FDA, you will update your labeling to reflect the additional testing. Such labeling updates will be made in consultation with, and require concurrence of, DMD/OHT7-OIR/OPEQ/CDRH.
- S. You must update the software application with the agreed upon FDA proposed edits within 1 month of the date of this letter (unless otherwise agreed to with DMD/OHT7-OIR/OPEQ/CDRH) and notify FDA that the agreed upon updates have been made prior to the distribution of your software application.
- T. You must evaluate the clinical performance of your product in children 13 years of age and younger in an FDA agreed upon post authorization clinical evaluation study within 6 months of the date of this letter (unless otherwise agreed to with DMD/OHT7-OIR/OPEQ/CDRH). After submission to and concurrence with the data by FDA, you must update the authorized labeling to reflect the additional testing. Such labeling updates will be made in consultation with, and require concurrence of, DMD/OHT7-OIR/OPEQ/CDRH.
- U. You must complete the agreed upon real-time stability study for your product and notify DMD/OHT7-OIR/OPEQ/CDRH of the testing results as they become available until completion of the study. After submission of the study data, and review and concurrence with the data by FDA, you must update your product labeling accordingly. Such labeling updates must be made in consultation with, and require concurrence of, DMD/OHT7-OIR/OPEQ/CDRH.
- V. You must evaluate the clinical performance of your product to support the serial screening claim in an FDA agreed upon post authorization clinical evaluation study within 6 months of the date of this letter (unless otherwise agreed to with DMD/OHT7-OIR/OPEQ/CDRH). After submission to and concurrence with the data by FDA, you must update the authorized labeling to reflect the additional testing. Such labeling updates will be made in consultation with, and require concurrence of, DMD/OHT7-OIR/OPEQ/CDRH.

⁸ Traceability refers to tracing analytical sensitivity/reactivity back to an FDA-recommended reference material.



Quick Aid
WHEN EVERY SECOND COUNTS

CELLTRION OTC RAPID ANTIGEN

Page 9 – Ally Danta, Celltrion USA, Inc.

- W. You must submit your product for any FDA-recommended independent evaluation to confirm the performance characteristics of your test, if requested by FDA. After submission to and concurrence with the data by FDA, you will update your labeling to reflect the additional testing. Such labeling updates will be made in consultation with, and require concurrence of, DMD/OHT7-OIR/OPEQ/CDRH.

Conditions Related to Printed Materials, Advertising and Promotion

- X. All descriptive printed matter, advertising, and promotional materials relating to the use of your product shall be consistent with the authorized labeling, as well as the terms set forth in this EUA and meet the requirements set forth in section 502(a), (q)(1), and (r) of the Act, as applicable, and FDA implementing regulations.
- Y. No descriptive printed matter, advertising, or promotional materials relating to the use of your product may represent or suggest that this test is safe or effective for the detection of SARS-CoV-2.
- Z. All descriptive printed matter, advertising, and promotional materials relating to the use of your product shall clearly and conspicuously state that:
- This product has not been FDA cleared or approved; but has been authorized by FDA under an EUA;
 - This product has been authorized only for the detection of proteins from SARS- CoV-2, not for any other viruses or pathogens; and,
 - The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

The emergency use of your product as described in this letter of authorization must comply with the conditions and all other terms of this authorization.

V. Duration of Authorization

This EUA will be effective until the declaration that circumstances exist justifying the authorization of the emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 is terminated under Section 564(b)(2) of the Act or the EUA is revoked under Section 564(g) of the Act.



Quick Aid
WHEN EVERY SECOND COUNTS

CELLTRION OTC RAPID ANTIGEN

Page 10 – Ally Danta, Celltrion USA, Inc.

Sincerely,

Jacqueline A. O'Shaughnessy, Ph.D.
Acting Chief Scientist
Food and Drug Administration

Enclosure